



IMAGE-PRO[®] EXPRESS

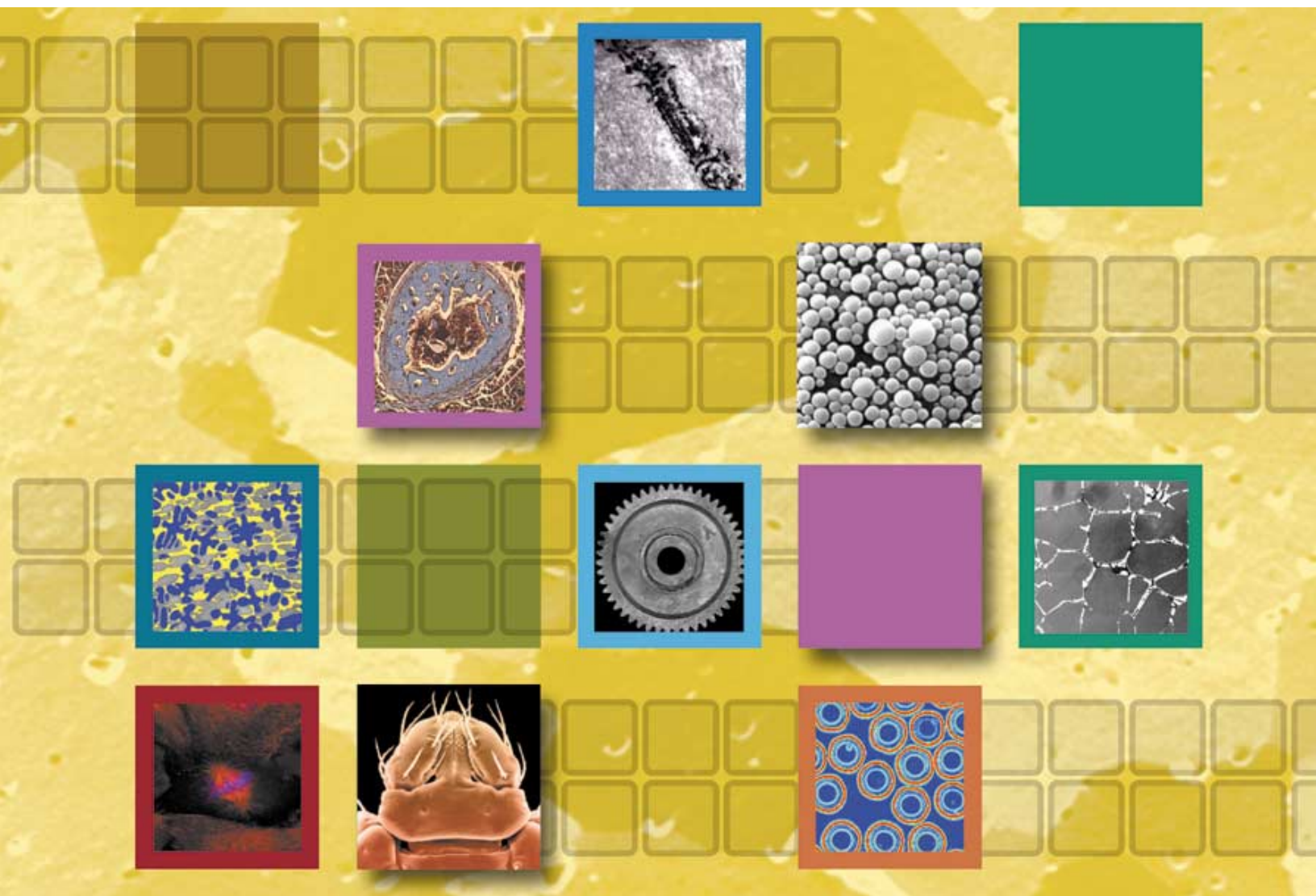


THE COST-EFFECTIVE IMAGE CAPTURE AND ENHANCEMENT SOLUTION FROM THE IMAGE-PRO FAMILY

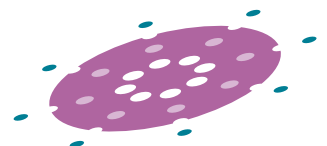

MediaCybernetics[®]
From Images to Answers[®]

Image-Pro[®] EXPRESS

Create in-focus composite images from a series of partially in-focus images with Extended Depth of Field (EDF)

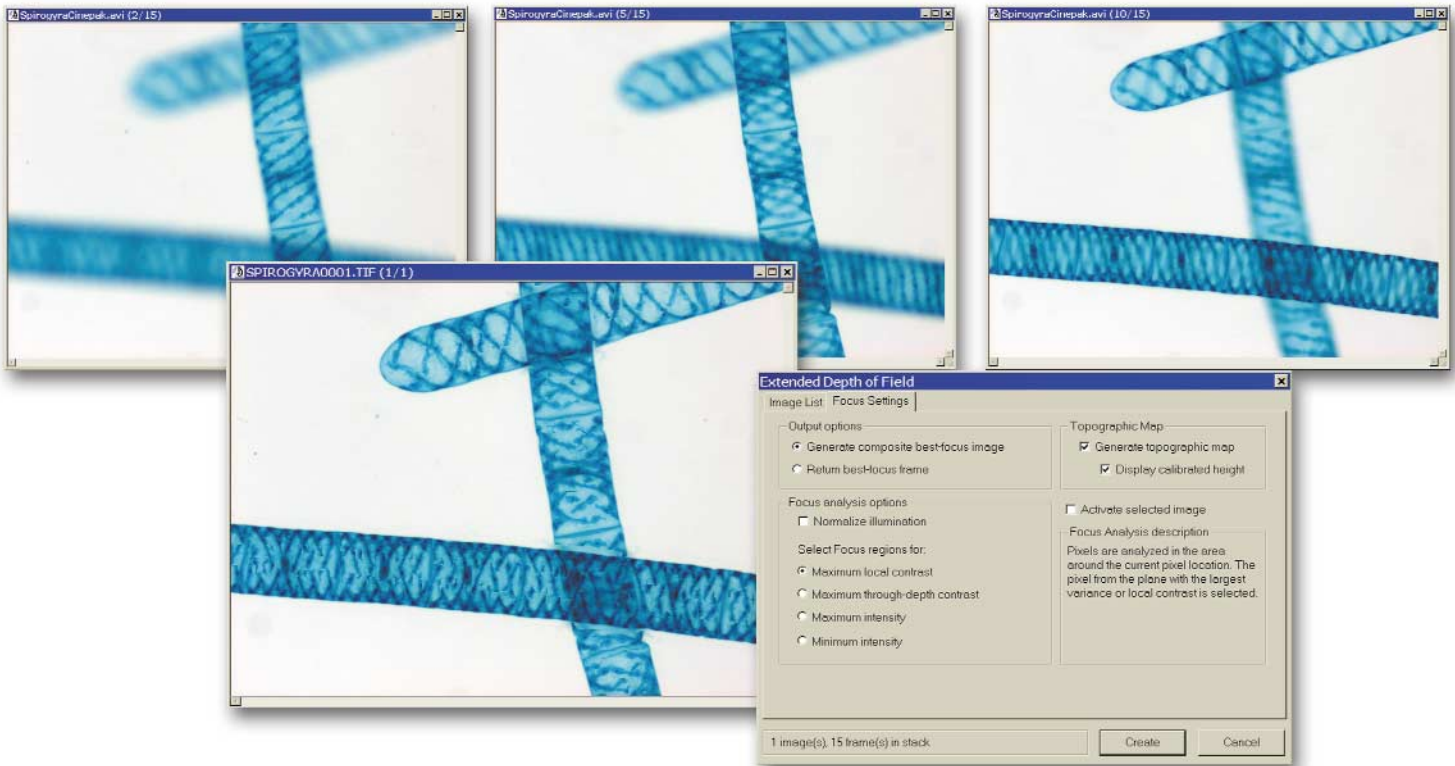

Experience the Speed and Quality of Image-Pro Express...

From Images to Answers

Image-Pro Express is the cost-effective image capture and enhancement software perfect for basic imaging or image capture stations. As the first step in the Image-Pro software series, Image-Pro Express comes equipped with the basic features needed to capture and enhance images for scientific, medical, and industrial research.

Numerous Capture Options

Image-Pro Express is a perfect addition to your imaging station. It allows you to quickly capture and enhance images with most popular digital cameras, video cameras, scanners, microscopes, and image databases. Images captured or saved at an Image-Pro Express station can also be analyzed with Image-Pro Discovery or Image-Pro Plus.

Easy Image Enhancement

Image-Pro Express allows you to quickly process images right at your capture station to guarantee high quality, accuracy, and reproducibility. Tools such as Extended Depth of Field (EDF) allow you to prepare images for detailed analysis. All enhancements made with Image-Pro Express are compatible with Image-Pro Discovery and Image-Pro Plus.

Upgrade for Additional Features

As your imaging requirements grow, you may find that you will need the power of Image-Pro Discovery or Image-Pro Plus. These

next steps in the Image-Pro Family offer advanced imaging tools and the opportunity to become part of an imaging community. You can upgrade to either of these products easily and at an affordable price.

The Image-Pro Family

The Image-Pro software series reflects the maturity of over 20 years of development, evolution, and user feedback. We offer three unique Image-Pro packages to accommodate a variety of imaging needs.

■ Image-Pro[®] Express

Perfect for basic imaging or for image capture stations, Image-Pro Express offers quick capture and enhancement tools.


■ Image-Pro[®] Discovery

More advanced than Image-Pro Express, Image-Pro Discovery offers added measurement and analysis capabilities.

■ Image-Pro[®] Plus

The ultimate imaging solution, Image-Pro Plus includes all of the features of Image-Pro Discovery along with added tools and the ability to write customized macros and plug-ins.

Apply Crosshair grids for precise placement of images during capture

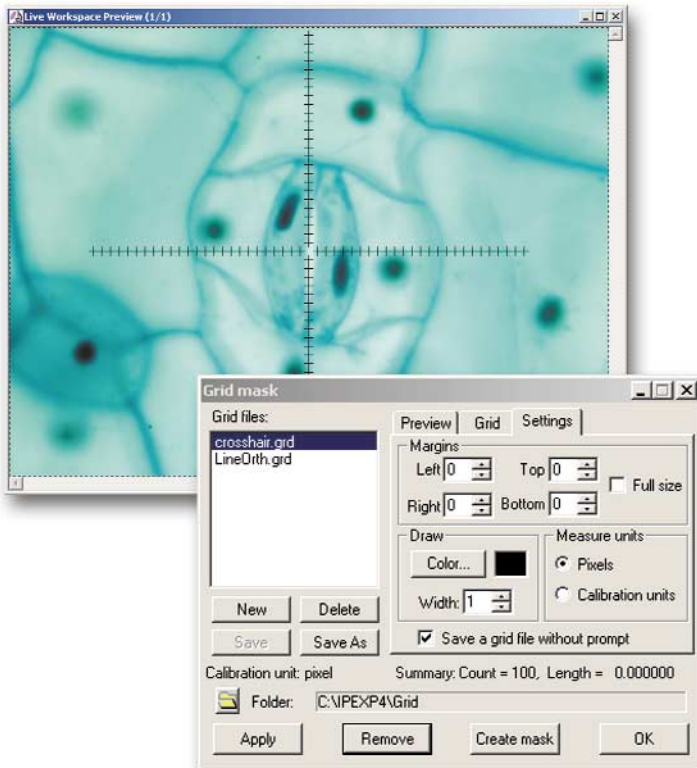

Measure images before you acquire them with Live Image Preview

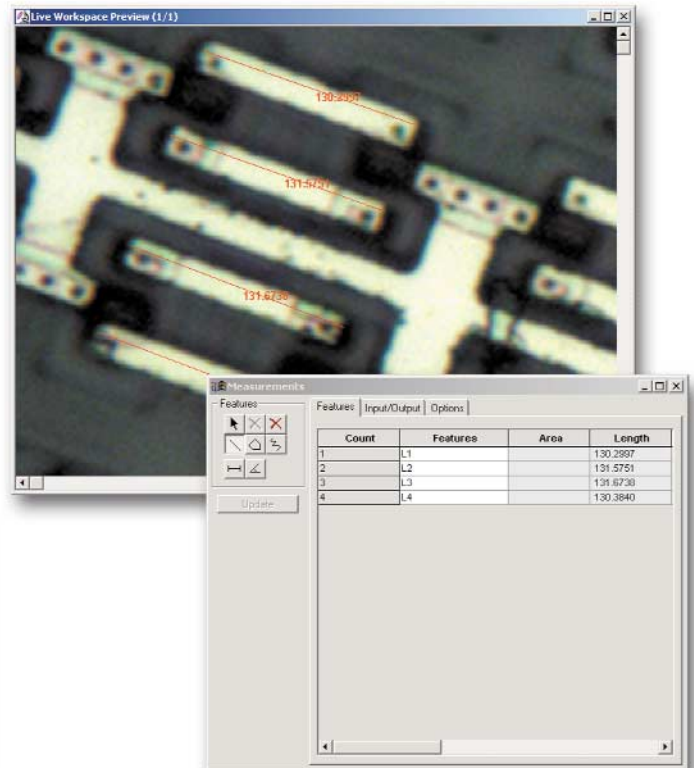


Image-Pro Express Features...

Acquire

- Acquire, create, and play back a sequence (stack) of images
- Acquire images from cameras, scanners, or electronically stored files
- Capture multiple images at user-specified time intervals

Enhance

- Edit image intensity level with Display Range
- Control image Brightness, Contrast, and Gamma with Contrast Enhancement
- Enhance with Best Fit, Linear, Bell, Logarithmic, and Exponential Equalization tools
- Invert, reset, and apply contrast to images
- Use automatic image inversion

Process

- Use Background Subtraction and Correction
- Work with Enhancement filters that include Low-Pass, Hi-Pass, Gaussian, Hi-Gaussian, Sharpen, Flatten, Median, Rank, and Local Equalization
- Create in-focus composite images from partially in-focus image stacks with EDF
- Apply filters to frames in a sequence
- Apply crosshair grids for precise placement of images during capture
- Extract and merge multiple color channels (including RGB, HSI, HSV, and YIQ) simultaneously with Color Channel tools
- Create, print, or save a collage of several images with Mosaic Image

Measure

- Calibrate intensity and spatial measurements
- Create and display spatial calibration markers
- Work with pre-defined spatial calibration units
- Outline objects automatically with Auto-Trace features
- Measure live images with Live Image Preview
- Export measurements to statistical and spreadsheet packages via DDE

Analyze

- Calculate straight line, circle, irregular line, or area histograms
- Display raw data or histograms
- Define and manage multiple areas of interest

Annotate

- Work with non-destructive image annotation
- Add text, simple graphics, and arrows
- Modify palette or colors


Generate Reports

- Create customized reports with images, data, and text
- Create and save report templates

Print and Publish

- Print frames from a sequence or montage of single images with Mosaic Printing feature
- Output high quality gray scale and color images to Windows® printers
- Experience complete control of image position on your page
- Print titles to make poster-sized prints
- Use Test Strip feature to select best output
- Use Screen Capture feature to capture an individual window, a user-defined area, or an entire screen to an image file

Acquire, create, and play back a sequence of images

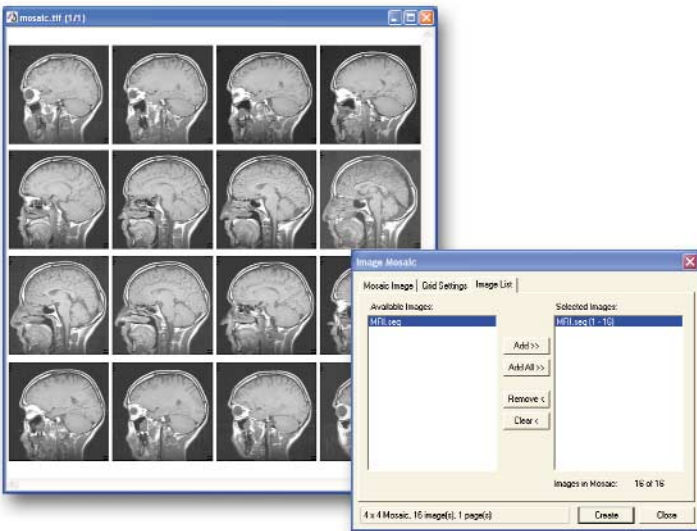

Create reports with images and data

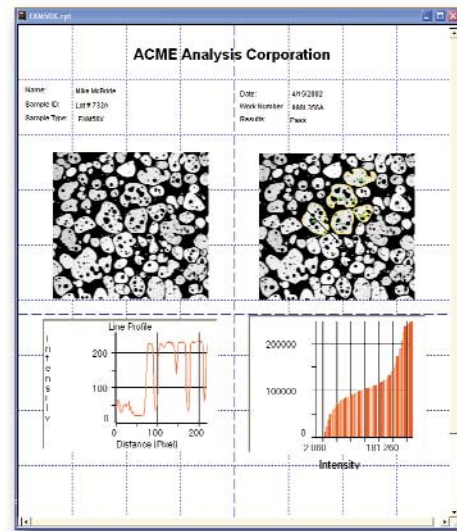


Image and Data File Format Support

- Read and write files in TIFF, JPEG, BMP, CUT, EPS, HDF, Photo CD, PCT, PCX, TGA, and FLB (Flat Binary)
- Read and write AVI files
- Convert files to and from the above formats
- Work with 8-, 12-, 16-, 24-, and 32-bit floating point images
- Compress images with JPEG, LZW, and RLE
- Batch convert files
- Output data files to ASCII, WKI, or XLS as input for spreadsheets
- Transfer images, data graphs, and data files via DDE or Clipboard
- Make best use of memory for images up to 64MB
- Store images on capture boards or in memory

System Requirements

Operating System

Operating System	Computer	RAM	Storage
Windows 98	133 Mhz or better	32 MB or more	8 MB
Windows ME	150 Mhz or better	32 MB or more	480 MB – 650 MB
Windows NT 4.0	200 Mhz or better	68 MB or more	650 MB
Windows 2000	133 Mhz or better	68 MB or more	650 MB
Windows XP	300 Mhz or better	128 MB or more	1.5 GB

- IBM Compatible PC
- A VGA (or better) resolution 256 color display device
- Windows compatible mouse, keyboard, and CD-ROM drive

OUR WORLDWIDE DEALER NETWORK IS THERE TO ASSIST YOU. FOR A DEALER IN YOUR AREA, CALL TODAY!

HEADQUARTERS

Media Cybernetics, Inc.
 8484 Georgia Avenue, Suite 200
 Silver Spring, MD 20910 USA
 Phone: +1.301.495.3305
 Fax: +1.301.495.5964
 Email: info@mediacy.com
 Web: www.mediacy.com

INTERNATIONAL OFFICES

South America:
 Phone: +55.11.4427.7803
 Fax: +55.11.4427.9527

Europe:
 Phone: +31.715.730.639
 Fax: +31.715.730.640

Asia Pacific:
 Phone: +65.6245.4965
 Fax: +65.6245.4967

